



Yes, you can.®

➤ Invacare® **Orion™**

**Stable, Trusty, Durable: The *Invacare® Orion™***

With a focus on safety and reliability, the **Orion™** is equipped with multiple features to ensure a safe and easy ride. Imagine being able to undertake long journeys in complete safety and comfort. Whether you enjoy a daily trip to the shops or love to travel further afield, the *Invacare® Orion™* gives you freedom of independence. Combining ergonomics with style, **Orion™** is the obvious choice for a relaxing, reliable drive.





### ➤ **Go further for longer**

Whatever the terrain, you'll enjoy a controlled and sturdy drive. Let the powerful motor, fully suspended chassis and big wheels take you smoothly over obstacles or uneven ground. The large batteries ensure plenty of capacity to enjoy excursions away from home. A unique lighting system includes an easily visible DUAL-bulb headlight, automatic indicators and an integrated heavy-duty brake light for safety. Lights operate in such a way that all bulbs function even if one light blows.

### ➤ **Enjoy relief and reliability**

Ergonomic steering on the **Orion™** is specifically designed to give relief to hands and arms during longer trips. User-friendly adjustment of the tiller is positioned opposite the handbrake and within easy reach. Together with the soft, adjustable swivel seat, this makes for an effortless, enjoyable drive. The **Orion™** is equipped with an electronic turning speed reduction feature ensuring safe, comfortable manoeuvres at all times.


### ➤ **Trustworthy and dependable**

Intelligent electronics offer a visual and audio warning when the batteries need charging.

A two-step disengaging lever prevents the scooter from free-wheeling if knocked out of drive mode. The **Orion™** has undergone extensive tests in the German independent TÜV test laboratory and complies with the Medical Device Standards for Europe.

### ➤ **Simplicity**

By only using top quality components, we have made the **Orion™** reliable and easy to service. Simply remove the top cover to access batteries, electronics, motor and cables. Electronics and the motor are protected from corrosion and water (splash protection) which safeguard the vital workings of the scooter against dirt and water preventing any unnecessary malfunction.

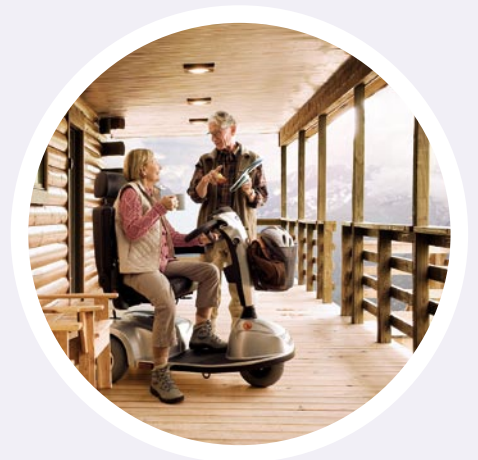

## Features



### Finite tiller movement

Users can easily adjust the tiller with a lever to suit their needs.



### Ergonomic steering

Prevents fatigue in arms and hands during long trips.



### Hand brake

Ensures immediate braking if required.



### Brake light

Heavy duty brake lighting warns people that the scooter is slowing down.



### Two-step disengaging lever

Prevents the scooter from accidentally free-wheeling.



### Built-in anti-splash guards

Protects electronics and transaxle from dirt and water.



### Automatic speed reduction

For added safety, the **Orion™** automatically reduces speed around bends.



### Suspension

A fully-suspended chassis provides a comfortable drive at all times.

## Options & Accessories



### Crutch & cane holder

Keeps your crutch or cane stable when on the move.



### Lockable front box

Secure box for personal valuables.



### Cup holder

Stores a soft drink.



### Heavy-duty storage cover

Keeps the scooter clean and dry.



### Off-board charging

Simply remove batteries to charge with the off-board charger kit.

## Technical Data

3 wheels	510 mm	470 mm	420 - 495 mm	510 mm	625 mm	1270 mm	103 kg	71 kg
4 wheels	510 mm	470 mm	420 - 495 mm	510 mm	625 mm	1290 mm	110 kg	78 kg
3 wheels	2 x 50 Ah	1125 mm	80 mm	90 mm	8°	42 km	10 km/h	10 km/h
4 wheels	2 x 50 Ah	1400 mm	80 mm	90 mm	12°	42 km	10 km/h	10 km/h

### Colour of shrouds



Moonstone Silver



Onyx Blue

### Upholstery



Black vinyl



10 km/h : 240 W / 600 W\*

\* Peak performance

All information quoted is believed to be correct at time of print. Invacare® reserves the right to alter product specifications without prior consultation.



INVACARE® Limited  
Pencoed Technology Park · Pencoed · Bridgend CF35 5HZ  
Tel: +44 1 656 776 222 · Fax: +44 1656 776 220  
E-mail: uk@invacare.com · [www.invacare.co.uk](http://www.invacare.co.uk)  
Sales Order E-mail: ordersuk@invacare.com

INVACARE® Ireland Ltd  
Unit 5 · Seatown Business Campus  
Seatown Road · SWORDS · County Dublin  
Tel. +353 1 8107084 · Fax +353 1 8107085  
Email: Ireland@invacare.com · [www.invacare.ie](http://www.invacare.ie)